




VariLase[®] Short

Laser Fiber and Procedure Kits

Another innovation from Vascular Solutions—Vari-Lase Short Kit

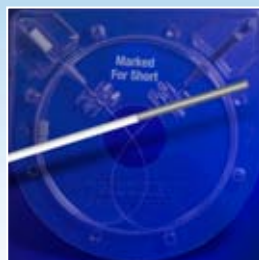

The Vari-Lase Short Kit is an endovenous laser kit designed specifically for the treatment of **short vein** segments

Designed to make access into these small vein segments easier, the Vari-Lase Short Kit features:

- Unique 0.018" guidewire with 1cm floppy tip for easy access and sheath support in short veins
- 10cm micro-introducer sheath for flexible, low-profile placement of Platinum Bright Tip fiber into short veins
- 21G needle in both standard 7cm and short 4cm length for ease of use with short veins
- Highly echogenic Platinum Bright Tip fiber specifically marked for use with micro-introducer sheath
- Fiber lock pre-positioned on Platinum Bright Tip fiber to conveniently secure the fiber to the micro-introducer sheath



VariLase® Short

Laser Fiber and Procedure Kits

Model	Description	Quantity
7161	Vari-Lase Platinum Bright Tip Short Kit <i>5F Micro-Introducer Set, 7cm needle and 600 micron micro-introducer marked Bright Tip laser fiber with pre-positioned fiber lock</i>	1
7282	10-Pack Short Kit Micro-Introducer Set <i>5F Micro-Introducer Sheath, 4cm smooth needle, 1cm floppy tip 0.018" 40cm nitinol mandrel</i>	1
7700	21G, 7cm smooth tip needle	1



The Vari-Lase Short kit is indicated for the treatment of varicose veins and varicosities associated with superficial reflux of the Great Saphenous Vein, and for the treatment of incompetence and reflux of superficial veins in the lower extremity. Please see the Vari-Lase Instructions for Use for a complete listing of indications, contraindications, warnings and precautions.

CAUTION: Federal law (U.S.A.) restricts these devices to sale by or on the order of a physician.

U.S. and International patents pending.

Vari-Lase is a registered trademark of Vascular Solutions, Inc.

Vari-Lase Bright Tip is a trademark of Vascular Solutions, Inc.

©2011 Vascular Solutions, Inc. All rights reserved.

ML1847 Rev. D 06/11



Vascular Solutions, Inc.

6464 Sycamore Court

Minneapolis, Minnesota 55369 USA

USA Customer Service: 888.240.6001

www.vasc.com